Pistolet do masažu- Massage gun-Masážní pistole-
Massagepistole

PDM3123

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
Item description

- 1) Battery removal button
- 2) Battery power indicator
- 3) Battery switch
- 4) Charging port
- 5) 6 kinds of massage heads
- 6) 99; 80; 60; 40; 20- % Power indicator
- 7) Gear indicator. For 01 gear, pulse frequency is 1200 rpm, if you add one gear the frequency will increase 100 rpm, the highest grade is 20, and the top speed is 3200 rpm.
- 8) Increase/ decrease buttons
- 9) Battery:(Voltage 24V)

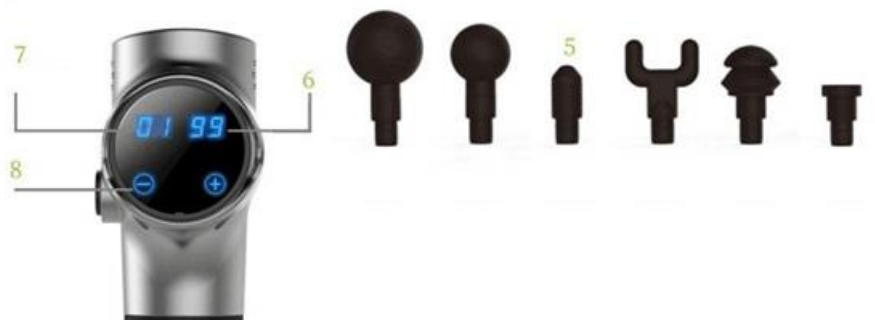Instructions for use:

1. To install the battery, please align the battery bulge to removal button, then insert the battery. To remove, please hold the battery grip and press battery removal button, then pull it out.
2. When you set the battery switch to "ON" position, the battery power indicator lights up.
3. According to the needs of different parts, select the appropriate massage head to insert into the round hole in front of the massager
4. Click the + button at the bottom of the screen-and the massager starts running ,the number on the top left of the screen represents the gear position, each + click button, increase the gear position by one gear,01 speed is 1200 rpm, every additional gear increased speed by 100 rpm, the highest gear is 20 gears and the speed is 3200 rpm. Once you click - button, the machine minus one gear. In any gear position (01-20), long press + or - button for five seconds will stop running.
5. After the massage is complete, turn off the battery switch to the "OFF" position and the screen will turn off.
6. In order to protect the machine, the massager will stop running automatically if it is running continuously for more than 10 minutes, or if the pressing force is too big, please turn it on again if you need to keep using it.

Battery use instruction

When the battery is exhausted, it needs 1.5-2.5 hours to be fully charged.

The battery working time is determined by the gear position. The battery power-on indicator lights up when the battery is charging, it will turn from red to green color.

For your safety, you should never turn on the machine when the battery is charging.

Safety instruction

To reduce the risk of electric shock, fire and personal injury, this product must be used in accordance with the following instructions.

1. Only allowed for adults use.
2. This device cannot directly massage the surface of the skin to avoid injury to the skin, the massage area should be worn after wearing dry and clean clothes or a clean and dry towel. You need to gently press and move, the massage time of the same part can not exceed 60s.
3. Use high speed and high pressure only for soft tissue parts of the body without causing pain and discomfort, do not use high speed and high pressure on the head and any bone parts of the body.
4. According to the massage needs of different parts, please choose the right massage head.
5. Frequent use of high-speed high pressure to apply pressure on the same part will cause abrasions, if you feel painful and uncomfortable, please stop using immediately.
6. Finger, hair or any part of the body, please keep away from the telescopic rod of the device and the heat dissipation hole of the motor to prevent from pinching injury.
7. Do not place anything in the heat sink of the motor.
8. Keep away from water
9. Do not discard and misuse the device
10. Use only original charger.
11. Carefully check your device and battery before use.
12. Do not make any changes in the device
13. Do not use the device and charge it in same time.

Maintenance

1. Wipe the surface of the device with a slightly damp towel and dry with a soft cloth.
2. When storing or traveling for a long time, separate the device from the battery and place it in the device's original case.

Warning:

1. Do not stay in the water, keep the machine away from liquids and heat sources, to ensure the motor cooling holes are free of dust and debris.
2. Do not remove the screws by yourself.
3. Turn off the device after one hour of operation and use it after rest for half an hour.
4. Remove the charger from the battery after charging, and do not use the device during charging. 5. The use of this device for the injured part of the body is prohibited without consulting the doctor .

Speed level:	1 ÷ 20
Vibration frequency (linear move):	10mm
Screen type:	LED
length (mm):	178
width (mm):	70
height (mm):	235
Motor power	20W
Output speed/ frequency	1200-3200 rpm
ZPower suply	230V
Battery:	24V
Power type	accumulator
Charging time	2,5-4h
Battery capacity	1500mAh
Working time after full charge	180min
Qty of massage head	6
Noise level	50dB
Included in set	Case
Material:	ABS; Silicon;
GW/NW	1,9/0,98KG

Labelling the product with a symbol of a crossed waste container informs us that the worn-out electrical and electronic equipment cannot be put into the same container as other waste. According to the WEEE Directive on the management of electrical and electronic waste, separate recycling methods shall be used for this type of equipment.

A user who wants to dispose of such product is obliged to transfer it to the collection centre for worn-out electrical and electronic equipment, thanks to which they support its re-usage, recycling or recovery and thusly protect natural environment. For this purpose, one should contact the equipment seller or representatives of local authorities. Hazardous components contained in electronic equipment may have

long-lasting harmful impact on natural environment as well as adverse effect on human health

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Date of sale:

GUARANTEE TERMS

1. The Seller provides guarantee on behalf of the Guarantor within the territory of the Republic of Poland for the period of 24 months from the date of sale.
In the case of sales and use product for the commercial purposes warranty period is equal 12 MSC from sales data
2. The Guarantee will be recognised by the shop or service centre after the client provides:
 - clearly and correctly filled-in guarantee card with the sale stamp and the seller's signature - valid purchase confirmation for the equipment including the date of sale / bill / - complained product.
3. Possible defects and damages revealed during the guarantee period shall be remedied free of charge within not more than 21 days from the delivery of the product to the shop or service.
4. Should it be necessary to import some components, the guarantee term can be extended by the period required to import such part, however such period shall not be longer than 40 days.
5. The guarantee does not cover:
 - mechanical damages and subsequent defects,
 - damages and defects resulting from incorrect storage and usage of the equipment against its purpose,- incorrect assembly and maintenance,
 - damages and consumption of such elements as: cables, straps, rubber parts, pedals, sponge grips, wheels, bearings, etc.
6. The guarantee expires in the event of:
 - expiry of the validity term,
 - self-repairs,
 - failure to follow the rules of correct operation.
7. No duplicates of the guarantee cards shall be issued.
8. The product handed over for repair should be complete and clean. If some parts are missing, the service is entitled to refuse to accept the product for repair. The service may refuse to accept a dirty product or clean in on the client's expense after the client's written permit.
9. The guarantee does not cover activities related to assembly, maintenance, which, according to the instruction, the user is obliged to perform themselves.
10. The Guarantor also informs that they offer post-guarantee service.
11. The product should be provided in its original packaging and should be secured for shipment.

The guarantee for sold goods is not exclusive, does not limit or suspend entitlements of the buyer resulting from the non-compliance of the product with the agreement.

Notes on the course of repairs

Item	Date of notification	Date of provision	Course of repairs	Signature of the recipient (shop, owner)

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